

EXAMION®

X-Ray Systems · Digital Imaging · Service

X-R Portable

Portable X-ray system



Customized solutions
for medical imaging

Why choose a portable X-ray system from EXAMION?

Your everyday practice routine has you often – or always – going to your patients' stables, or the animals are too big to be examined at your practice? Then you need X-ray equipment which comes along with you and can be used flexibly. The portable X-ray solutions by EXAMION can help you with this task.

- Energy-saving LED light for collimator
- Convenient remote control via hand switch
- Flat touch keyboard
- Digital LED display
- Display reversal for improved readability from various angles.
- Robust and dirt-repellent housing
- Programmable speed buttons for frequently used kV and mAs values
- Control of source image distance via dual laser



Battery-powered devices

EXAMION battery-powered devices are the best choice for mobile veterinary practice. They combine a battery-powered device with an X-DR portable solution, consisting of the flat panel and notebook. Get perfectly matched components and a wireless mobile X-ray for a complete examination.



X-R Portable Air Ray

The X-R Portable AirRay device is the lightest among the X-ray generators. The powerful battery allows up to 120 exposures per battery charge and can even be connected by cable.

- Weight: 5.3 kg
- Performance: max. 90 kV, max. 20 mA, 1.6 kW bei 80 kV

X-R Portable 20 BT

The X-R Portable 20 BT is the most compact device. The very low weight is an essential feature to reduce the physical workload during pre-purchase examinations.

- Weight: 7.2 kg
- Performance: max. 90 kV, max. 20 mA, 1.6 kW @ 80 kV



X-R Portable 20 BT Plus

The X-R Portable 20 BT plus is a perfect choice for continuous operation. It has a larger housing and effectively dissipates the heat that is generated during imaging.

- Weight: 11.2 kg
- Performance: max. 100 kV, max. 20 mA, 1.6 kW @ 80 kV

X-R Portable 40 BT

The X-R Portable 40 BT is the perfect choice if power is a crucial feature.

- Weight: 14.2 kg
- Performance: max. 100 kV, max. 35 mA, 2.4 kW @ 80 kV



Wired devices

The hybrid equipment from EXAMION can be used for stationary and mobile applications. The hybrid equipment with cables is suitable for users who mainly work at their practice but would occasionally like to use the equipment for a mobile application. Equipment with cables also has advantages when it is rarely used and the user does not wish to deal with checking and recharging the battery.

X-R Portable 20 Plus

The X-R Portable 20 Plus is the perfect choice for the mobile horse specialist.

- Weight: 9.8 kg
- Performance: max. 100 kV, max. 20 mA, 1.6 kW @ 80 kV



X-R Portable 40

The X-R Portable 40 is the perfect choice for both stationary and mobile use. The right choice for veterinarians covering both small animals and mobile horse applications.

- Weight: 12.8 kg
- Performance: max. 100 kV, max. 35 mA, 2.4 kW @ 80 kV



X-More: Accessory products by EXAMION

Portable solutions – everything inside

The components of compact X-ray solutions are ideally matched and provide an optimal supplement to our portable X-ray devices.

- Detector
- Laptop or Ultrabook
- Carrying case or tote bag
- Image acquisition using our X-AQS



X-ray tripods

An X-ray tripod is perfectly suited for convenient working when taking X-rays – particularly with mobile equipment. It provides relief for the user and makes X-ray examinations more efficient.



X-R MORE Stand Stat-x portable

- Very quick and convenient assembly without tools
- Flexible and height adjustable
- Very lightweight (only 15 kg), but can be loaded with up to 25 kg
- When collapsed:
30 cm x 30 cm x 142 cm
- Available in three models:
 1. Collapsible mobile stand for portable x-ray equipment
 2. For CR cassettes or detectors
 3. Combination of options 1 and 2

EXAMION GmbH

EXAMION is a global provider of digital radiography solutions for human and veterinary medicine. We have been operating according to our, one-stop expert' principle for more than 30 years, assisting our customers with competent advice in planning, installation and after-sales services.



We currently support more than 7,000 systems, which are primarily stationary installations in hospitals and medical practices, partially also mobile applications. Our proprietary EXAMION X-AQS software for system control, image processing, diagnosis and archiving ensures a standardised control and efficient workflow.


In addition to radiography and digital image processing, our product range includes mammography and radiotherapy.

Our more than 150 specialists work as a team for our customers throughout Germany, and four European locations in Belgium, Poland, Portugal and Switzerland.

Through a network of certified sales and service partners, EXAMION provides customers around the world with its latest state-of-the-art X-ray solutions.

EXAMION is characterised by its extensive market expertise, cutting-edge technology and its vast number of project references in the private and public sector.

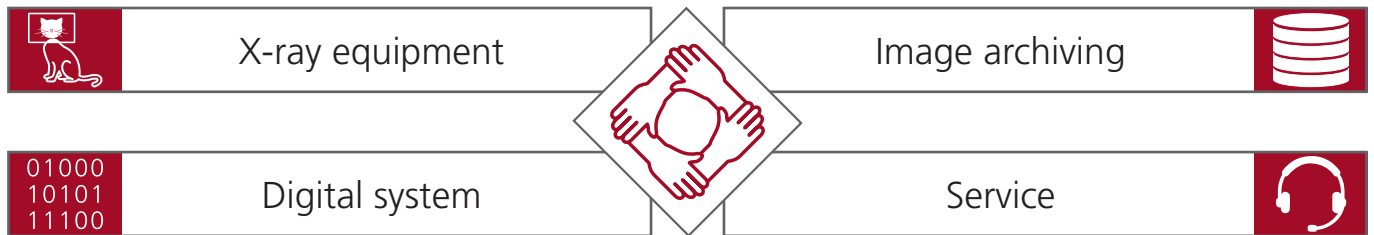


Benefit from years of experience



Contact head office

EXAMION GmbH
Erich-Herion-Str. 37
70736 Fellbach
Germany

www.examion.com

 0711 120 002 - 0
 0711 120 002 - 22
 sales@examion.com



Certified partner of EXAMION